



Comparative study of displacement resistance of four zirconia cements

Estudio comparativo de la resistencia al desplazamiento de cuatro cementos en zirconia

Enrique Ríos Szalay,* Alfredo Garcilazo Gómez,^{§,||} Jorge Guerrero Ibarra,*
Ivonne Meade Romero,* Karla Miguelena Muro*

ABSTRACT

Objective: To compare displacement resistance of four cementing agents. **Material and methods:** An experimental, cross-sectioned prospective research was conducted to assess four cementing agents. Three agents were resinous, self-adhesive, dual polymerization cements containing MDP (10-metacryloxydecyl dihydrogen phosphate), and the remaining was a conventional glass ionomer cement. In the experiment, 40 samples of zirconia partially stabilized with yttrium were prepared. All samples were treated following their specific manufacturer's instructions. Samples were prepared, they were then stored at 100% humidity in a temperature chamber at 37 °C for 24 hours; after this, samples were subjected to shearing detachment mechanical tests at a 1 mm per minute speed in a universal machine for mechanical testing. **Results:** Glass ionomer samples failed before being taken to the universal testing machine. Remaining three cements did not show statistically significant differences. **Conclusions:** Adhesion capacity of glass ionomer to zirconia is nil or extremely low. Likewise, resinous cements containing MDP in their formula, either in their bonding agent or in the cement formulation itself, are presently the best alternative to increase adhesion to a zirconia structure.

Key words: Cementing agents for zirconia.

Palabras clave: Agentes cementantes para zirconia.

RESUMEN

Objetivo: Comparar la resistencia al desplazamiento de cuatro agentes cementantes. **Material y métodos:** Se realizó una investigación prospectiva, transversal y experimental en la que se evaluaron cuatro agentes cementantes, tres de ellos resinosos autoadhesivos de polimerización dual y con contenido de MDP (10-metacrilóxido decil dihidrógeno fosfato) y un ionómero de vidrio convencional. Se realizaron 40 muestras de zirconia parcialmente estabilizada con itrio, se dividieron en cuatro grupos, cada uno de ellos fue tratado de acuerdo con las indicaciones del fabricante del cemento a estudiar, se realizaron las muestras, se almacenaron en humedad al 100% en una cámara a una temperatura de 37 °C durante 24 horas para después ser sometidas a pruebas mecánicas de desprendimiento por cizallamiento a una velocidad de 1 mm por minuto en la máquina universal de pruebas mecánicas. **Resultados:** La muestras de ionómero de vidrio fracasaron antes de ser llevadas a la maquina universal, entre los otros tres cementos no existe diferencia estadísticamente significativa. **Conclusiones:** La capacidad de adhesión de ionómero de vidrio a la zirconia es nula o muy baja. Igualmente los cementos resinosos que contengan en su fórmula MDP, ya sea en su agente de acoplamiento o en la fórmula misma de los cementos, son en la actualidad la mejor alternativa para incrementar la adhesión a una superficie de zirconia.

INTRODUCTION

Zirconium oxide partially stabilized with yttrium (Y-TZP) better known as zirconia, has constituted a great success in the field of biomaterial research. Since the decade of the '70s, use of zirconia in dentistry was evidenced through studies proposing its use as a coating for implants.¹ Nevertheless, it was only in the '90s when there were first reports of its use in implants.² In 1991, there were reports of zirconia use in orthodontic brackets.³ Use of zirconia in the field of restorative dentistry began during the middle of this decade when it was used for manufacture of

* Graduate and Research School. National School of Dentistry, National Autonomous University of Mexico.

§ Health and Biological Science Division, Healthcare Department, Metropolitan Autonomous University, Xochimilco campus.

|| Division of Professional studies. National School of Dentistry, National Autonomous University of Mexico.

Received: January 2017.

Accepted: May 2017.

© 2017 Universidad Nacional Autónoma de México, [Facultad de Odontología]. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

This article can be read in its full version in the following page:
<http://www.medigraphic.com/facultadodontologiaunam>

intra-root posts and crowns manufactured with CAD/CAM as well as appliances for rehabilitation of dental implants and fixed partial prostheses.⁴⁻⁶ To the present date, zirconia treatments, due to their high values of fracture resistance, have become ideal candidates to manufacture ceramic prostheses in areas of high mechanical compromise.

The main attribute of Y-TZP (Yttria tetragonal zirconia polycrystalline) was described by Garvie in 1975, when he described the **resistance to transformation** phenomenon, in which, partially stabilized zirconia in tetragonal phase, in the presence of a high stress area such as the extreme of a crack, suffers phase change in that area, passing to crystallizing that area in the monoclinic phase. This change involves an approximately 5% volume increase of the zirconia particle, able to seal the crack. Thus, healing of the area is ultimately achieved arresting crack increase (*Figure 1*).⁷

Y-TZP is a fracture-resistant material with excellent mechanical properties, it is considered to be biotolerable, and provides flexural strength of more than 900-1200 MPa, these are values two to three times higher than maximum mastication forces (200 to 400 N in anterior teeth and up to 600 N in posterior teeth). This flexural strength is higher than that exhibited by any other previously developed ceramic materials for dental use.⁸ It also exhibits a yield strength higher than almost all metallic alloys used in dentistry, its elasticity module (205 GPa) is somewhat lower than that exhibited by stainless steel (210 GPa) and similar to that of titanium alloys (Ti6Al4V);⁹ it presents thermal conductivity lower than alumina (zirconium 2.5 W 7 Mk versus alumina 30 W7mk at 37

°C),⁸ therefore, probability of triggering hypersensitivity in the case of sudden thermal changes is decreased.

It is a highly biotolerable material¹⁰ with low radioactivity, with radio-opacity similar to that of metals,¹¹ allowing thus excellent radiographic contrast.

Nevertheless, zirconia is not devoid of problems, among them we can count spontaneous degradation (related to hydro-thermal transformation) and stress derived from manufacturing process.¹² With respect to an ideal cementing agent, even though many research projects have been conducted, to this date, there are no strong results to help us determine which cementing system can be more suitable or more effective, therefore, protocols with resinous cements as well as glass ionomer protocols are recommended.^{13,14}

Zirconium is an acid-resistant ceramic material, differing from vitreous porcelains, it does not react to acid etching, moreover, it is quite unstable when subjected to thermal and mechanical changes.¹⁵ Traditional protocols of acid etching with hydrofluoric acid and silanization used to adhere other ceramic structures to dental structure are not applicable to zirconia, since there is absence of vitreous matrix and its nature is relatively inert; this renders it a low reactivity surface.^{16,17} Development of selective acid etching methods, sanding or infiltration have been attempted in order to prime zirconium surfaces to chemically or micromechanically adhere to dental structure with the use of resinous cements, targeting improvement of their mechanical properties without generating stress on the structure which might cause fractures and thus lead to failure.^{14,18} Nevertheless, to the present date, there are no studies to support effectiveness and durability of new protocols proposed for roughness generation (sanding, three-fold mechanical/chemical treatment, porcelain pearls, plasma spray) and thus chemically activate the zirconia surface (silanization, acrylization, silicon tetrachloride vaporization, MPD silanes and cements).¹⁹

Presently the most widely used technique to cement zirconia restorations would be use of sanding with aluminum oxide micro-spheres (50-110 μ , 2 to 3 pressure bars, 3 to 4 cm distance) along with cementing agents which contain phosphate monomers (MDP)^{8,19} are perhaps the technique more frequently used to cement zirconia restorations. It has been shown that cements containing monomer 10 metacryloxydecyl dihydrogen phosphate (MDP) exhibit particular affinity to metallic oxides such a zirconium dioxide, alumina and metal. MDP is a relatively hydrophobic monomer, due to its 10 carbons chain; it contains a hydrophil phosphate terminal which chemically adheres to zirconium oxide, and a polymerizable methacrylate terminal which adheres to resin.²⁰

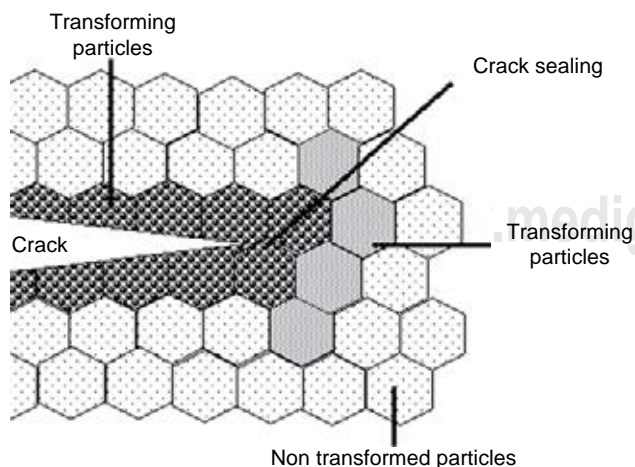


Figure 1. Representation of stress-induced transformation resistance process.

It is important to point out that tooth preparation is paramount among multiple factors leading to success of fixed restorations, as a factor of great influence in the retention degree of the restoration irrespectively of used material. Preparation of an abutment with convergence angles oscillating from 4° and 15°^{21,22} is widely found in the descriptions of clinical procedures of several follow up studies of zirconia restorations. It has been shown that decreasing the preparation's convergence degree to 10° exponentially increases retention degree,²³ irrespectively of used cement.

METHOD

Resistance to displacement of four cementing agents was compared, out of them, three agents were self-adhesive, dual polymerization cements with MDP content, and the remaining one was a conventional, glass ionomer cement. All cements were commonly used at the final stages of a zirconia rehabilitation. Cements were compared in order to determine which one would represent the best option for the aforementioned process.

The following cements were studied (*Figure 2*):

RelyX™ Ultimate with Single Bond® Universal, 3M™ ESPE™.

Multilink® Automix with Monobond® Plus, Ivoclar Vivadent.

PANA VIA™ SA Cement Automix, Kuraray Noritake Dental Inc.

Ketac Cem, 3M™ ESPE™.

Forty 7 x 7 mm square samples of Zirconia Lava™ Plus, 3M™ ESPE™ were obtained. They were sintered at 1,450 °C for 8 hours, according to manufacturer's instructions in an oven program S1P1600, Ivoclar Vivadent.



Figure 2. Cements used in the study.

Zirconia samples were soaked with PMM (methyl polymethacrylate), in 25 mm diameter polypropylene rings, using a different color for each study group (*Figure 3*). All samples were sanded with 50µ aluminum oxide spheres; as part of the cleansing process; samples were taken to a Branson 2510 ultrasound appliance for one minute. A Teflon shaper with a 4 mm orifice was used; with a press, cements to be tested were placed on the zirconia. In the first group, Multilink® Automix (Ivoclar Vivadent) was placed, after having applied bonding agent Monobond® Plus, of the same commercial brand.

Group 2 was previously treated with Single Bond® Universal (3M™ ESPE) as previous bonding agent, same procedure as previous group was undertaken with cement RelyX Ultimate 3M™ ESPE™.

Group three was achieved with PANAVIA™ SA Cement Automix, Kuraray Noritake Dental Inc.

Group four was undertaken with Ketac Cem (3M™ ESPE™).

All cements were handled according to manufacturer's instructions. Light-polymerizing cements were cured with an Ultradent Products Inc lamp at 600 mW/cm² power, measured with a Demetron brand radiometer and following timing established by manufacturer of each cementing agent. Specimens were then stored at 100% humidity in a chamber at 37 °C for 24 hours. In each sample area was calculated with the formula $\pi \times r^2$, area of all specimens was obtained in this manner. Mechanical test of dislodgment by shearing was conducted at a 1 mm per minute speed, in order to observe adhesion strength, a universal machine for mechanical tests Instron® model 5567 USA was used (*Figure 4*).

RESULT ANALYSIS

Once obtained, results were subjected to variance analysis test of one factor (ANOVA). Informative package SigmaStal® was used.

Ketac Cem (3M™ ESPE™) glass ionomer samples were discarded, since during manufacturing, all 10 samples failed (*Figure 5*).



Figure 3. Samples of all three resin cements.

According to standard deviation, the study was reliable for the following cements: Multilink® Automix Ivoclar Vivadent. Relyx™ Ultimate 3M™ ESPE. PANA VIA™ SA Cement Automix (Kuraray Noritake Dental Inc). Average values of retention forces of 7.223 MPa, 11.024 MPa and 12.256 MPa, respectively were reported (Table I).

ANOVA test with $p = 0.170$ to compare displacement resistance of the three cements revealed that there was no statistically significant difference with value $p = 0.070$. *Post hoc* tests were additionally performed, and likewise, it could be observed there were no statistically significant differences among groups, even though there was a 5 MPa difference range among cements of greater and lesser retention.

Null hypothesis was accepted based on results obtained in the present research project which had proposed that «there is no difference in resistance to displacement among cementing agents».

DISCUSSION

Based on data obtained from the tests, it was decided to exclude glass ionomer from the statistical analysis. Ernest et al (2005),¹⁶ Marchan et al (2005)²⁴ and Uo et al (2006)²⁵ reported that conventional ionomer showed very low retention values when used with zirconia. When comparing 5.8 MPa PANA VIA with Ketac Cem, Shahin and Kern reported results that showed that PANA VIA 21 exhibited best behavior as adhesion agent in the retention of zirconia structures without surface treatment and reported approximately 2.8 MPa.²⁶

Results of the present research confirm that null hypothesis; results obtained were similar to those

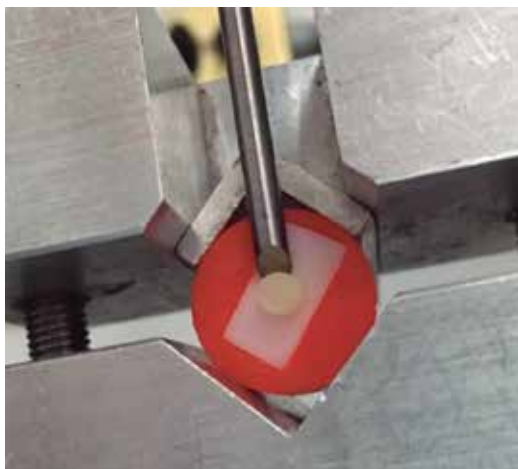


Figure 4. Shear test in the universal mechanical test appliance Instron® model 5567 USA.

of Palacios et al (2006)¹⁵ who tested three cements (PANA VIA F 2.0, Kuraray; RelyX Luting, 3M ESPE and RelyX Unicem, 3M ESPE), when they reported bonding to zirconia with no statistically significant differences.

Cements used in the present work are resinous cementing agents containing MDP in their formula. Kern and Wenger, in 1998²⁷ were the first to report long term adhesion strength of MDP-containing resinous cements; this was later confirmed in numerous studies.²⁸⁻³¹

Bonding agents have been developed in recent years. They have been introduced to improve bonding strength of ceramic to zirconium. When using RelyX® Elite and Multilink® Automix, the manufacturer recommends, in an alternative manner, to use systems with imprinters or bonding agent Single



Figure 5. Ketac Cem, 3M™ ESPE™ glass ionomer.

Table I. Mean and standard deviation of cements' retention strength.

Cement	N	Failures	Media	Standard deviation	SEM
Multilink	10	0	7.223	3.221	1.019
Ultimate	10	0	11.024	4.365	1.380
PANA VIA	10	0	12.256	6.402	2.024

SEM = Standard error of means.

Bond® Universal, 3M™ ESPE™ and Monobond® Plus, Ivoclar Vivadent respectively. In both cases, imprinters contain dihydrogen phosphate 10-methacryloxydecyl (MPD) and silanium. Numerous studies such as those of Amaral et al, Ozcan et al and Yoshida et al support the use of imprinters with MDP on zirconium oxide surfaces, since obtained results have shown that phosphate monomers are securing chemical agents in order to improve bonding to zirconia.^{20,30,31}

It has been concluded that roughness and activation of zirconia are important in order to achieve bonding of resin to the restoration. Many researchers (Bopna, Kern, Blatz among others), used abrasion techniques with aluminum oxide particles sanding on the restorations surface so as to increase surface energy, adhesion area and humectability.³²

McLaughlin (1984) Cortes (2003 and 2010) reported that either way and regardless of preparation to the zirconia surface, it is important to know that there will be no «integration» or «fusion» of restorations to the dental structure, as would be the case when using restorations with vitreous phase, treated with hydrofluoric acid and later silanized.³³⁻³⁵

CONCLUSIONS

Within limitation of the present study and taking into account it was conducted both in vivo and in vitro, it is possible to conclude the following:

- Glass ionomer adhesion capacity to zirconia is nil or extremely low.
- Presently, the best alternative to increase adhesion to a zirconia surface are resinous cements containing MDP in their formulation, either in their bonding agent or in the cement formula itself.
- To this date no statistically significant differences have been found among the three studied cements. PANAVIA offers a simplified placement technique, that is to say, it does not require bonding agents, thus decreasing probability of failure during cement manipulation.
- No studies have been found that support satisfactory cement adhesion to zirconia or in a similar manner to that achieved by vitreous phased ceramics.
- It is important to mention that in addition to technique or treatment applied to zirconia surfaces, long term success of the restoration is related to application of basic principles in the design preparation.
- We recommend conducting long term in vivo research projects in order to observe behavior of cementing agents used in zirconia oxide restorations.

REFERENCES

1. Cranin AN, Schnitman PA, Rabkin SM, Onesto EJ. Alumina and zirconia coated vitallium oral endosteal implants in beagles. *J Biomed Mater Res.* 1975; 9 (4): 257-262.
2. Minamizato T. Slip-cast zirconia dental roots with tunnels drilled by laser process. *J Prosthet Dent.* 1990; 63 (6): 677-684.
3. Springate SD, Winchester LJ. An evaluation of zirconium oxide brackets: a preliminary laboratory and clinical report. *Br J Orthod.* 1991; 18 (3): 203-209.
4. Meyenberg KH, Lüthy H, Schärer P. Zirconia posts: a new all-ceramic concept for nonvital abutment teeth. *J Esthet Dent.* 1995; 7 (2): 73-80.
5. Luthardt RG, Sandkuhl O, Reitz B. Zirconia-TZP and alumina--advanced technologies for the manufacturing of single crowns. *Eur J Prosthodont Restor Dent.* 1999; 7 (4): 113-119.
6. Yildirim M, Edelhoff D, Hanisch O, Spiekermann H. Ceramic abutments--a new era in achieving optimal esthetics in implant dentistry. *Int J Periodontics Restorative Dent.* 2000; 20 (1): 81-91.
7. McMeeking RM, Evans AG. Mechanics of transformation-toughening in brittle materials. *J Am Ceram Soc.* 1982; 65 (5): 242-246.
8. Montagna F, Barbesi M. *Cerámicas, zirconio y CAD/CAM.* México: Amolca; 2013. p. 151-169.
9. Piconi C, Rimondini L, Cerroni L. *El zirconio en odontología.* México: Amolca; 2011. p. 18-39.
10. Piconi C, Maccauro G. Zirconia as a ceramic biomaterial. *Biomaterials.* 1999; 20 (1): 1-25.
11. Porstendörfer J, Reineking A, Willert HC. Radiation risk estimation based on activity measurements of zirconium oxide implants. *J Biomed Mater Res.* 1996; 32 (4): 663-667.
12. Chevalier J. What future for zirconia as a biomaterial? *Biomaterials.* 2006; 27 (4): 535-543.
13. Manicone PF, Rossi Iommetti P, Raffaelli L. An overview of zirconia ceramics: basic properties and clinical applications. *J Dent.* 2007; 35 (11): 819-826.
14. Blatz MB, Sadan A, Kern M. Resin-ceramic bonding: a review of the literature. *J Prosthet Dent.* 2003; 89 (3): 268-274.
15. Palacios RP, Johnson GH, Phillips KM, Raigrodski AJ. Retention of zirconium oxide ceramic crowns with three types of cement. *J Prosthet Dent.* 2006; 96 (2): 104-114.
16. Ernst CP, Cohnen U, Stender E, Willershausen B. *In vitro* retentive strength of zirconium oxide ceramic crowns using different luting agents. *J Prosthet Dent.* 2005; 93 (6): 551-558.
17. Ozcan M, Vallittu PK. Effect of surface conditioning methods on the bond strength of luting cement to ceramics. *Dent Mater.* 2003; 19 (8): 725-731.
18. Aboushelib MN, Kleverlaan CJ, Feilzer AJ. Selective infiltration-etching technique for a strong and durable bond of resin cements to zirconia-based materials. *J Prosthet Dent.* 2007; 98 (5): 379-388.
19. Echeverri DM, Garzón H. Cementación de estructuras para prótesis parcial fija en zirconia. *Rev Fac Odontol Univ Antioq.* 2013; 24 (2): 321-335.
20. Yoshida K, Tsuo Y, Atsuta M. Bonding of dual-cured resin cement to zirconia ceramic using phosphate acid ester monomer and zirconate coupler. *J Biomed Mater Res B Appl Biomater.* 2006; 77 (1): 28-33.
21. Vult von Steyern P, Carlson P, Nilner K. All-ceramic fixed partial dentures designed according to the DC-Zirkon technique. A 2-year clinical study. *J Oral Rehabil.* 2005; 32 (3): 180-187.
22. Tinschert J, Schulze KA, Natt G, Latzke P, Heussen N, Spiekermann H. Clinical behavior of zirconia-based fixed partial dentures made of DC-Zirkon: 3-year results. *Int J Prosthodont.* 2008; 21 (3): 217-222.

23. Jørgensen KD. The relationship between retention and convergence angle in cemented veneer crowns. *Acta Odontol Scand.* 1955; 13 (1): 35-40.
24. Marchan S, Coldero L, Whiting R, Barclay S. *In vitro* evaluation of the retention of zirconia-based ceramic posts luted with glass ionomer and resin cements. *Braz Dent J.* 2005; 16 (3): 213-217.
25. Uo M, Sjögren G, Sundh A, Goto M, Watari F, Bergman M. Effect of surface condition of dental zirconia ceramic (Denzir) on bonding. *Dent Mater J.* 2006; 25 (3): 626-631.
26. Shahin R, Kern M. Effect of air-abrasion on the retention of zirconia ceramic crowns luted with different cements before and after artificial aging. *Dent Mater.* 2010; 26 (9): 922-928.
27. Kern M, Wegner SM. Bonding to zirconia ceramic: adhesion methods and their durability. *Dent Mater.* 1998; 14 (1): 64-71.
28. Amaral R, Ozcan M, Valandro LF, Balducci I, Bottino MA. Effect of conditioning methods on the microtensile bond strength of phosphate monomer-based cement on zirconia ceramic in dry and aged conditions. *J Biomed Mater Res B Appl Biomater.* 2008; 85 (1): 1-9.
29. Ozcan M, Nijhuis H, Valandro LF. Effect of various surface conditioning methods on the adhesion of dual-cure resin cement with MDP functional monomer to zirconia after thermal aging. *Dent Mater J.* 2008; 27 (1): 99-104.
30. Tanaka R, Fujishima A, Shibata Y, Manabe A, Miyazaki T. Cooperation of phosphate monomer and silica modification on zirconia. *J Dent Res.* 2008; 87 (7): 666-670.
31. Wegner SM, Kern M. Long-term resin bond strength to zirconia ceramic. *J Adhes Dent.* 2000; 2 (2): 139-147.
32. Magne P, Paranhos MP, Burnett LH Jr. New zirconia primer improves bond strength of resin-based cements. *Dent Mater.* 2010; 26 (4): 345-352.
33. McLaughlin G. Porcelain fused to tooth--a new esthetic and reconstructive modality. *Compend Contin Educ Dent.* 1984; 5 (5): 430-435.
34. Corts JP. *Carillas o frentes estéticos y sus variantes.* In: Lanata EJ y col. Atlas de operatoria dental. Buenos Aires: Alfaomega; 2008. p. 251.
35. Corts JP. *Restauraciones indirectas adheridas anteriores.* In: Henostroza G. Adhesión en Odont Rest de ALODYB. Madrid, España: Ripano SA; 2010. p. 34.

Mailing address:
Enrique Ríos Szalay
E-mail: enriosza@unam.mx